

Ames Archive Manager for the Ames Study Manager

The Ames Archive Manager offers extensive searching and filtering capabilities, helping you find specific studies for archiving quickly and easily.

The Ames Archive Manager belongs firmly within our software suite for regulatory Ames testing. The program has been developed in accordance with the OECD series on principles of GLP and compliance monitoring number 15: *Establishment and Control of Archives that Operate in Compliance with the Principles of GLP*.

Reconstruct archived studies for reporting, auditing and verification

Lock down study data whilst retaining flexibility to view or report on original data

Easily browse, search and filter all of your Ames studies, past and present

Developed with reference to OECD guideline principles of GLP and compliance monitoring number 15

GLP and FDA 21 CFR Part 11 compliance with audit trails, user permissions and database security

Comprehensively tested and validated in accordance with our Quality Management System which is approved to ISO 9001:2008 incorporating ISO 90003:2004



A safe and secure archiving utility for those seeking GLP compliance


The Ames Archive Manager is aimed at those wishing to archive their assay data in line with the OECD guideline on the 'Establishment and Control of Archives that Operate in Compliance with the Principles of GLP'. The program works in conjunction with the Ames Study Manager and its associated modules to provide the complete solution for regulatory Ames testing.

The program offers extensive searching and filtering capabilities, helping you find specific studies for archiving quickly and easily. Once archived, a study cannot be modified in any way, yet retains accessibility for report reconstruction, auditing and verification.

Archive studies with a single click without the need to export study data

The Ames Archive Manager allows users to select individual studies or a range of studies to be archived in just one simple action.

In order to maintain the integrity and accessibility of your data, the Ames Archive Manager works by simply marking selected studies as archived, and so does not remove any data from your original Ames database. This not only helps ensure the integrity of your data by avoiding any potentially dangerous export and purge operations, but also means you can still view the study and reconstruct reports for auditing and verification.

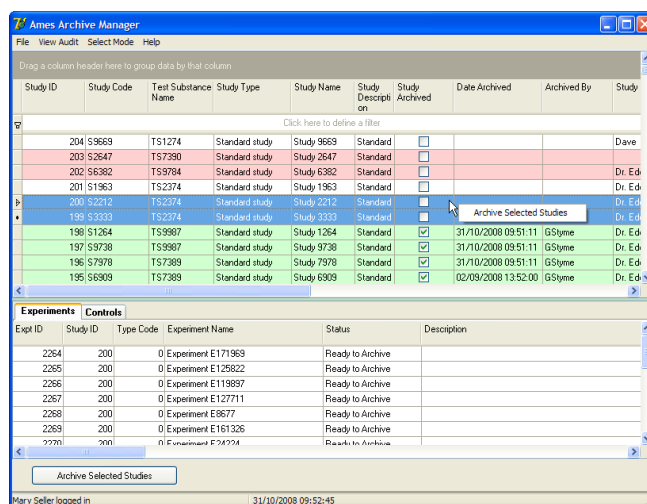

Study ID	Study Code	Test Substance Name	Study Type	Study Name	Study Description	Study Archived	Date Archived	Archived By	Study
204	S9668	TS1274	Standard study	Study 9668	Standard	<input type="checkbox"/>			Dave
203	S2847	TS7390	Standard study	Study 2647	Standard	<input type="checkbox"/>			
202	S6382	TS9784	Standard study	Study 6382	Standard	<input type="checkbox"/>			Dr. Ed
201	S1963	TS2374	Standard study	Study 1963	Standard	<input type="checkbox"/>			Dr. Ed
200	S2212	TS2274	Standard study	Study 2212	Standard	<input type="checkbox"/>			Dr. Ed
199	S3550	TS2274	Standard study	Study 3550	Standard	<input type="checkbox"/>			Dr. Ed
198	S1264	TS9987	Standard study	Study 1264	Standard	<input checked="" type="checkbox"/>	31/10/2008 09:51:11	GSyme	Dr. Ed
197	S9739	TS9987	Standard study	Study 9739	Standard	<input checked="" type="checkbox"/>	31/10/2008 09:51:11	GSyme	Dr. Ed
196	S7978	TS7389	Standard study	Study 7978	Standard	<input checked="" type="checkbox"/>	31/10/2008 09:51:11	GSyme	Dr. Ed
195	S6309	TS7389	Standard study	Study 6309	Standard	<input checked="" type="checkbox"/>	02/09/2008 13:52:00	GSyme	Dr. Ed

Expt ID	Study ID	Type Code	Experiment Name	Status	Description
2264	200	0	Experiment E171969	Ready to Archive	
2265	200	0	Experiment E125822	Ready to Archive	
2266	200	0	Experiment E119897	Ready to Archive	
2267	200	0	Experiment E127711	Ready to Archive	
2268	200	0	Experiment E8677	Ready to Archive	
2269	200	0	Experiment E161326	Ready to Archive	
2270	200	0	Experiment E26774	Ready to Archive	

Auditing and user administration helping you meet regulatory requirements

The Ames Archive Manager has been designed and developed with reference to a wide range of regulatory guidelines and when used in conjunction with the System Access Manager is fully compliant with Good Laboratory Practice and the FDA 21 CFR Part 11 final rule on electronic records and electronic signature.

The program has been tested and validated in accordance with our Quality Management System which is approved to ISO 9001:2008 incorporating ISO 90003. Comprehensive audit trails record every action taken within the program, including both the archiving and unarchiving (where allowed) of study data. These audit events can be viewed not only from within the Ames Archive Manager itself, but also from the main Ames Study Manager main event log.

